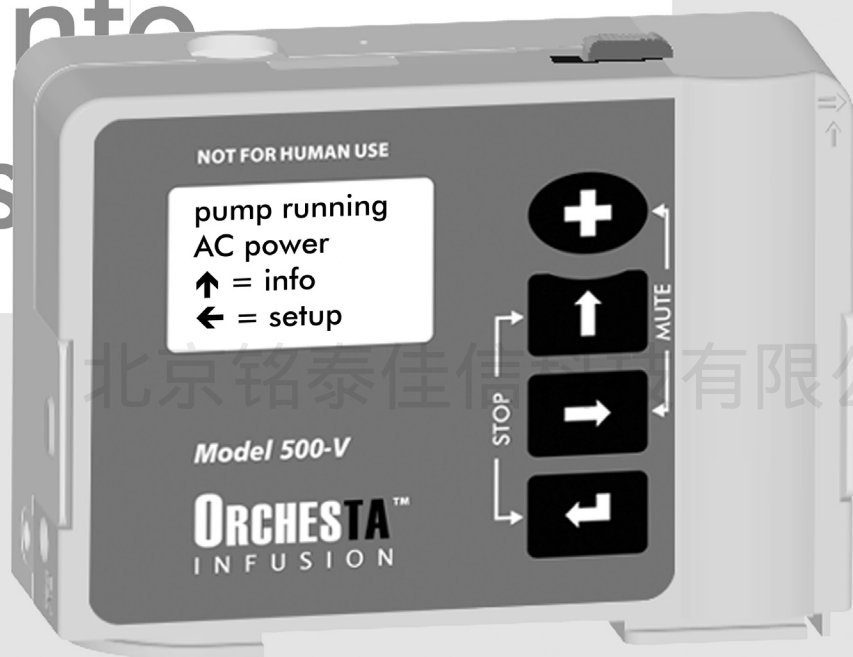


Instruction manual



for serial numbers after LV-25000

Infusion pump

ORCHESTA™ Model 500-V

CE

NOT FOR HUMAN USE



ORCHESTA is a trademark of Instech Laboratories, Inc.

manufactured by:

exclusive lab animal distribution by:



Mühlenstraße 17
D - 24229 Dänischenhagen b. Kiel
Phone +49 (0) 4349 91 54 0
Fax +49 (0) 4349 91 54 10
info@vennermedical.de

INSTECH

Instech Laboratories, Inc.
5209 Militia Hill Road
Plymouth Meeting, PA 19462 USA
www.instechlabs.com
Tel (800) 443-4227
Tel +1-610-941-0132
Email support@instechlabs.com

Contents

Introduction	4
Design and functional overview	8
Preparation for infusion	15
Power supply	18
Programming	24
Brief instructions for daily care	36
Information menu	42
Alarm and warning functions	43
Cleaning and safety check	51
Dosing accuracy confirmation	52
Compliance and classification	53
Technical data and accessories	54
Appendix	56
Index	58

Introduction

Designated Use

ORCHESTA™ Model 500-V is a portable, volumetric micro-infusion pump for the administration of fluids for laboratory animal use only, not for human use. It can be used for multiple routes of administration.

The main accessories for the infusion pump are: *PEGA® PUMP HEAD V*, *PEGA® BAG*, *PEGA® LINE*, *PEGA® SPIKE* and *PEGA® BOX*.

PEGA® BAGs serve as containers for fluids and are connected when full to the inlet side of the *PEGA® PUMP HEAD V* in the pump via a Luer-Lock connection.

The *PEGA® BOX* is provided for the *PEGA® BAG* in sizes of 50 ml, 100 ml and 150 ml.

PEGA® BOX with the *PEGA® BAG* is slid onto the side of the pump thereby creating a stable connection between the two.

PEGA® LINE connects the pump to the animal. *PEGA® SPIKE* is used with infusion bags having a septum.

PEGA® PUMP HEAD V, *PEGA® BAG*, *PEGA® LINE* and *PEGA® SPIKE* are all disposables.

Notes on the use of the instruction manual

This manual is intended for technically instructed staff (**see p. 5, Items 2 and 3**).

International safety symbols are used in this instruction manual. Paying attention to these symbols will protect you from injury and avoid hazards. Please also note the key to the symbols on the back of this instruction manual.

Warning

Disregarding these instructions may result in serious damage to animal health or material damage.

Caution

Disregarding these instructions may result in material damage.

Note

Additional information, which is not related to safety.

General Safety Information

1. Read the manual before setting up the system and be sure to observe the safety instructions.
2. ORCHESTA™ Model 500-V may only be used with the supervision of personnel whose knowledge and judgement define the suitability of the pump for a particular infusion application. Disregarding this requirement will invalidate the warranty (**see p. 56**).
3. Using and changing the PEGA® PUMP HEAD V, PEGA® BAG and PEGA® catheter systems (e.g. PEGA® LINE) may only be carried out by qualified staff.
4. Do not operate the pump higher than 1 m above the animal's heart.
5. Carry out the prescribed safety checks regularly. If use of the pump is continued beyond the next safety check date, the warranty will be invalidated and all product liability of the manufacturer will lapse.
6. In order to prevent risks, use only original PEGA® consumable materials and accessories (**see p. 55**).
7. Check that the pump is working and the unit is in proper condition each time before using it. To do this, a self-test must be carried out when starting up (**see p. 22**) with priming the pump (**see p. 33**).
8. PEGA® PUMP HEAD V, PEGA® BAG and PEGA® catheter systems (e.g. PEGA® LINE) are disposable items in sterile packaging. Only use if the packaging is undamaged and unopened.
9. For optimal aseptic technique, consider replacing the PEGA® PUMP HEAD V, PEGA® BAG and PEGA® catheter systems (e.g. PEGA® LINE) at the same time.
10. Only connect filled PEGA® BAGs to the pump in the primed state.
11. Dispose of disposable items according to local regulations.
12. During mains/AC operation, infusion therapy may be interrupted due to mains/AC-borne interference. Such interference may occur if the network is overloaded by a lightning strike for example. In such a case, please follow the instructions in

the operating manual under the item "Start-up/Self-test" (**page 22/23**) to restart the infusion.

13. The field strength of stationary radio transmitters should be in compliance at all frequencies according to on-site tests. Interference is possible in the vicinity of equipment which carries the adjacent graphical symbol.



14. Do not use zinc carbon batteries as these may cause the unit to malfunction. Note the battery /

rechargeable battery options shown in **Table 1 on p. 19**.

15. Never expose ORCHESTA™ Model 500-V to X-rays (e.g. CT), strong magnetic fields (e.g. MRI) or ultrasound.
16. Defective units or units that are no longer required must not be disposed of in the household waste. Observe the specific disposal regulations for electrical equipment for your country.
17. Batteries and rechargeable batteries should be disposed of properly.

Additional information

- ◆ ORCHESTA™ Model 500-V has a splash-proof housing meaning that it will tolerate a small amount of water exposure.
- ◆ If a malfunction occurs, the pump will set off a visual and audible alarm, if enabled.
- ◆ The safety concept of the pump is designed in such a way that it

reverts to the safe state if a fault occurs:

If a fault occurs, infusion is stopped and a visual and audible alarm is signaled.

If the power supply is interrupted during dosing, there is only an audible alarm. The visual alarm is switched off.

北京铭泰佳信科技有限公司

Design and functional overview

Hardware

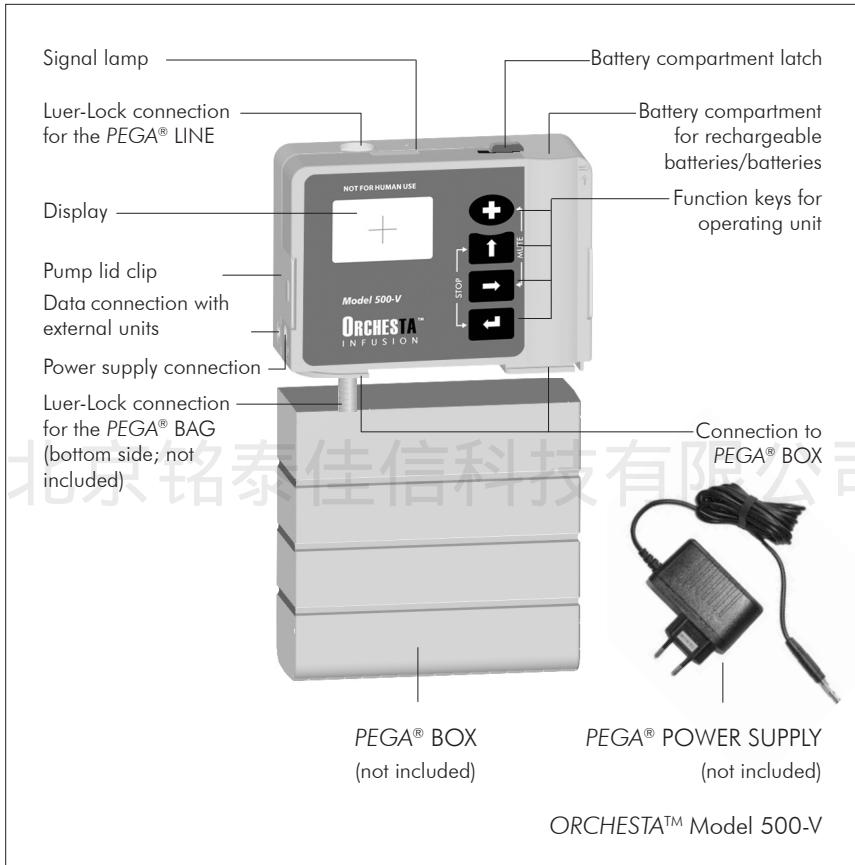


Figure A: View of pump, PEGA® BOX, and PEGA® BAG assembled

Principle of operation

The ORCHESTA™ Model 500-V infusion pump delivers fluids from a reservoir (e.g. PEGA® BAG) through a catheter

system (e.g. PEGA® LINE) to the animal. In the filled state, the reservoir is connected to the bottom of the pump by

means of a Luer-Lock connection, and the catheter system is connected to the PEGA® PUMP HEAD V by means of the Luer-Lock connection on the top of the pump (see p. 16ff). While pumping is in progress, a micro-piston drive, which is

driven by a motor in the housing, ensures that the liquid is correctly pumped. Rechargeable or ordinary disposable batteries are used to provide power to the pump. The PEGA® POWER SUPPLY can also be used.

Operating modes

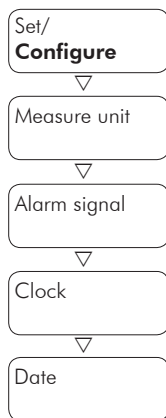
This manual will guide you through the continuous infusion mode only.

Operational Menus

The ORCHESTA™ Model 500-V has three operational menus. The settings menu 1 and 2 are secured by codes. The information menu is freely accessible. The information menu differs depending on the unit of measurement set in settings menu 1. If in settings menu 1 no key is pressed for more than 30 min-

utes or the power supply is interrupted, all parameters must be reentered.


If in the information menu no key is pressed for about 15 seconds, the screen backlight switches off in order to save power. It can be switched on again by pressing any key.

Settings menu 1**Code 1 - 9 - 2****Sub-menu A (Configure)****Selection:** Configure**Unit of measure:** ml/h, mg/h, $\mu\text{g/h}$,
mg/d or $\mu\text{g/d}$ **Alarm signals:** select between different alarm signals (internal warning tones cannot be changed)**Clock:** sets the current time**Date:** sets the current date

北京铭泰佳信科技有限公司

Settings menu 1

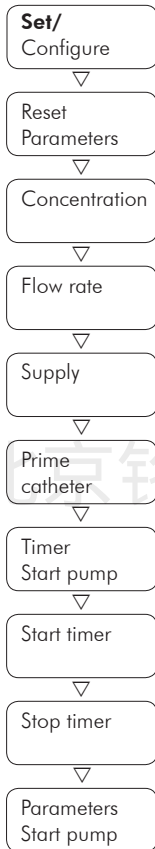Code 1 - 9 - 2



Sub-menu B (set)

Selection: Set

Reset: deletes all parameters set in this menu

Concentration: mg/ml, $\mu\text{g/ml}$ (is only shown if mg/h, mg/d, $\mu\text{g/h}$ or $\mu\text{g/d}$ is set as the unit of measurement in sub-menu A "Configure")

Flow rate: rate at which fluid is continuously administered

Supply: amount of fluid with which the PEGA® BAG is filled; 1 ml - 1,000 ml in steps of 1 ml

Prime catheter: filling the catheter system with the fluid; flow rate during priming: 15 ml/h


Start timer: set timer for scheduled pump start

Stop timer: set timer for scheduled pump stop

Start pump: terminates settings menu 1 and starts the pump with the parameters set

Settings menu 2

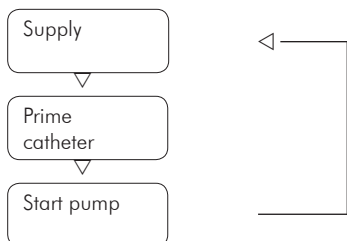Code 1 - 1 - 1




Supply: enter new amount of article after the PEGA® BAG has been changed

Prime catheter: filling the catheter system with the administration article; flow rate during priming: 15 ml/h

End: terminates settings menu 2, the pump returns to its previous state

Information menu

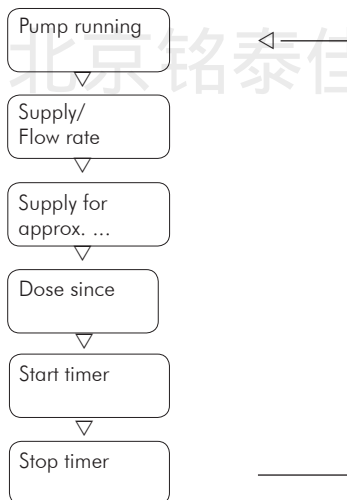No code



Unit of measure in ml/h

Supply/Flow rate: current fluid amount in the PEGA® BAG/flow rate set

Supply for approx.: range of fluid supply in the PEGA® BAG

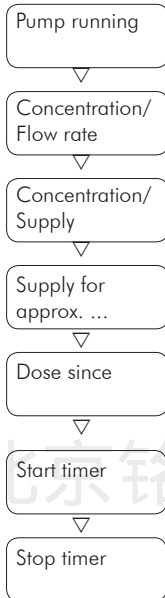
Dose since...: amount of fluid administered during the last 24 hours

Start timer: time for scheduled pump start (only if programmed)

Stop timer: time for scheduled pump stop (only if programmed)

Information menu

No code



Unit of measure in mg/h, mg/d, $\mu\text{g/h}$ or $\mu\text{g/d}$

Concentration/Flow rate: concentration of the fluid (mg/ml or $\mu\text{g/ml}$)/set flow rate

Concentration/supply: concentration of the fluid (mg/ml or $\mu\text{g/ml}$)/current amount of the fluid in the PEGA® BAG

Supply for approx.: range of fluid supply in the PEGA® BAG

Dose since....: amount of fluid administered in the last 24 hours

Start timer: time for scheduled pump start (only if programmed)

Stop timer: time for scheduled pump stop (only if programmed)

Function

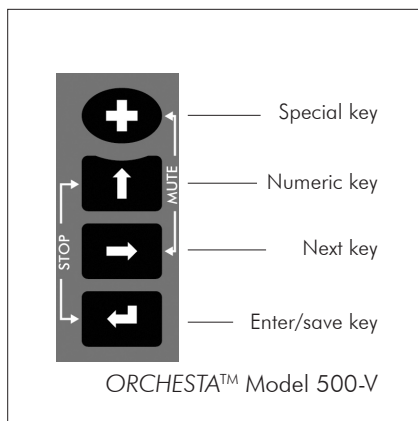


Figure B: Function keys

Special key: primes the catheter, re-starts an interrupted program

Numeric key: enters numerical values

Next key: displays the next message in the information or settings menu, selects the next numeral, prompts the information menu

Enter/save key: changes or saves/deletes the entry, prompts the window for entering the code

Stop: Pressing the **Numeric key** and the **Enter/Save key** simultaneously stops the pumping process.

Mute/Unmute: Pressing the **Special key** and the **Next key** simultaneously in the Information menu mutes/unmutes the audible signals.

Transport and storage

Transport

Appropriately cover and protect the ORCHESTA™ Model 500-V when placing in an animal jacket. When transporting the pump, use the plastic case with dedicated foam cutouts in which it was shipped originally.

Storage

If you intend not to use ORCHESTA™ Model 500-V for a long period of time, stop the pump, switch it off (**see p. 40**), remove the batteries or rechargeable batteries and store ORCHESTA™ Model 500-V in its own case. The data will remain stored until the next safety check.

Preparation for infusion

Open infusion pump

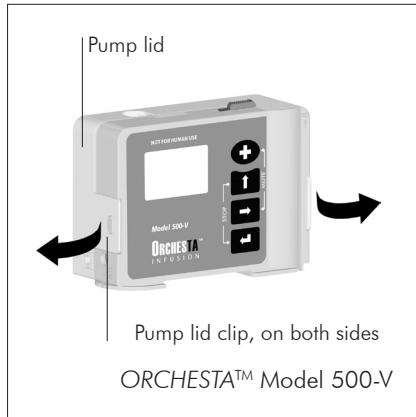


Figure C: Side view

Warning

The pump may only be opened for changing disposables. Further dismantling of the pump housing may lead to restrictions in the operation of the pump and will invalidate the warranty.

Procedure

1. Disconnect the PEGA® BOX if it is connected.
2. Place the ORCHESTA™ Model 500-V on a flat surface.
3. Slightly bend the two pump lid clips outwards, one after the other.
4. Lift the pump lid off away from you.



Figure D: Side view with pump lid open

Inserting the PEGA® PUMP HEAD V

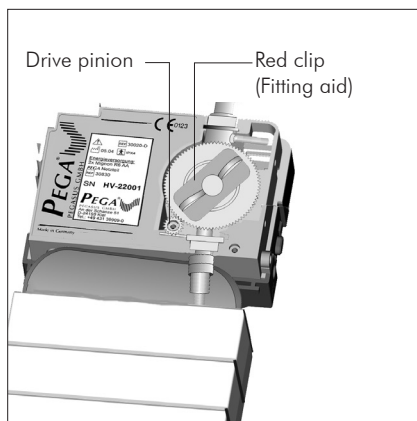


Figure E: PEGA® PUMP HEAD V - Fitting

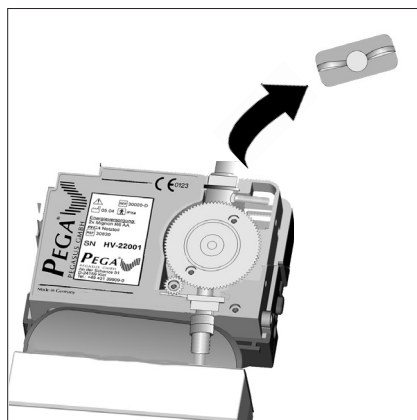


Figure F: PEGA® PUMP HEAD V – Removing the clip

Procedure

1. It is essential to follow the general safety instructions (**see page 5**).
2. Open the pump.
3. Open the sterile packaging of the PEGA® PUMP HEAD V.
4. Hold the PEGA® PUMP HEAD V by the red clip (fitting aid).
5. Insert the PEGA® PUMP HEAD V into the pump so that the drive pinion slides over the drive shaft at the bottom left (**Figure E**).
6. Make sure that the gears are not misaligned when inserting the PEGA® PUMP HEAD V!
7. Press with the fingers on the Luer connectors until the PEGA® PUMP HEAD V is firmly seated in the pump.
8. Carefully remove the red clip (**Figure F**).
9. Place the pump lid on the pump.
10. Close the pump. Make sure that the lid clips on the side engage fully into place.

Note

ORCHESTA™ Model 500-V is only ready for operation and will only operate properly if the clips are correctly engaged.

⚠ Caution

In order to avoid contamination, do not open the protective cap of the PEGA® PUMP HEAD V until the PEGA® BAG and PEGA® LINE are

connected. Do not touch the Luer connections after removing the protective cap.

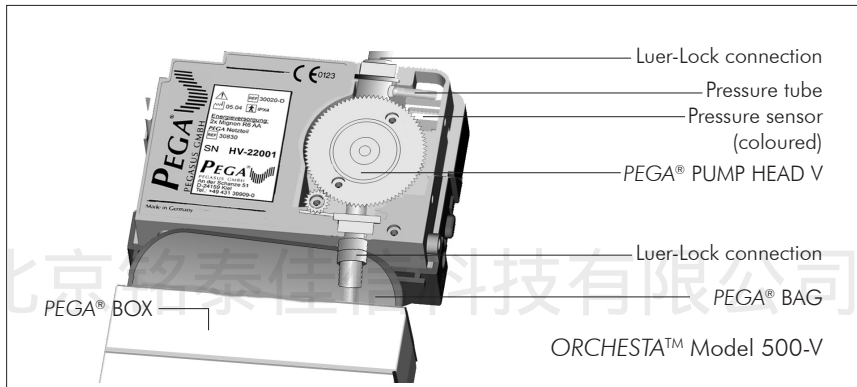


Figure G: Rear view without pump lid

PEGA® Disposables

To connect the PEGA® BAG and PEGA® catheter system (e.g. PEGA® LINE), please refer to the instructions enclosed

with the disposable materials. You will find the steps needed for the pump operation on **page 36ff.**

PEGA® BOX and soft cases

To accommodate the 50 ml, 100 ml and 150 ml PEGA® BAG connect the corresponding PEGA® BOX with the pump. Solomon Scientific offers stable pouches

for using the infusion pump without the PEGA® BOX. They have been specially designed for this application and ensure a safe delivery without the PEGA® BOX.

Power supply

To guarantee correct operation of the ORCHESTA™ Model 500-V, the pump should be operated using the batteries provided. However, it is also possible to use batteries (ordinary, disposable) or rechargeable batteries (Type R6, AA).

The PEGA® POWER SUPPLY is available for operating the pump from the mains/AC supply. Solomon Scientific also provides an external charger with rechargeable batteries.

The life of normal batteries and rechargeable batteries (Type R6, AA) depends on their capacity and the flow rate set. The resulting approximate values for battery life are shown in Table 1 (battery/storage battery life) on **page 19**.

Notes

- ◆ The pump is unable to detect rechargeable batteries or normal batteries which are inserted during operation. Therefore disconnect from mains/AC power before inserting and confirming new batteries or rechargeable batteries.

- ◆ If the batteries have been removed by mistake before the pump has been stopped or switched off, an audible alarm will sound which can only be halted by re-inserting the batteries/rechargeable batteries or by connecting to AC power.

Caution

- ◆ Use only the separately available PEGA® POWER SUPPLY for supplying the unit with AC power.
- ◆ Do not use zinc carbon batteries as this may affect reliable operation of the pump.
- ◆ Rechargeable batteries and normal batteries have different voltage levels. After inserting, you must tell the pump which type of cell has been used when requested.

This is necessary in order to guarantee correct operation of the pump and the internal monitoring functions and alarms.

Inserting rechargeable battery pack/batteries, connecting the power supply unit

Batteries	Capacity in mAh	at basal rate of 7.5 ml/h		at basal rate of 1 ml/h	
		Range	Life	Range	Life
Lithium	3000	200 ml	27 hrs.	150 ml	150 hrs.
Alkali-Manganese	2250	150 ml	20 hrs.	100 ml	100 hrs.
NiMH rechargeable pack PEGA® RECHARGEABLE	1650	110 ml	15 hrs.	80 ml	80 hrs.
NiCd rechargeable	1000	70 ml	9 hrs.	50 ml	50 hrs.

Table 1: Options of batteries/rechargeable batteries



Figure H: Battery compartment lid open

Fitting/changing rechargeable batteries and normal batteries

Procedure

1. Pull out the AC plug if necessary.
2. Place the pump on a flat surface with the battery compartment lid upwards.
3. Release the latch of the battery compartment lid, push the lid up and move it sideways.
4. Remove any empty or defective batteries or rechargeable batteries.
5. Insert two new 1.5 V batteries or rechargeable batteries (Type R6, AA). Pay attention to the polarity (label in the battery compartment).
6. Fold the lid back over the compartment, gently press it downward and close the latch.

Operating with power supply unit

Procedure

1. Connect the power supply unit (PEGA® POWER SUPPLY) to the mains/AC.
2. Plug the power supply unit into the socket on the infusion pump.
3. The unit is now working on AC power.

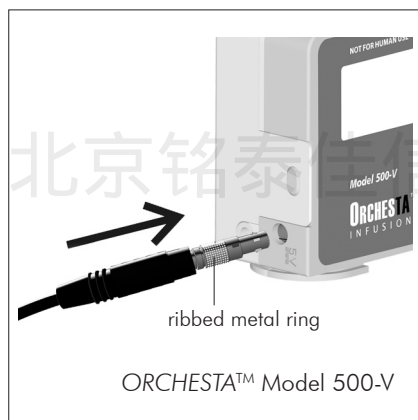


Figure 1: Socket for power supply

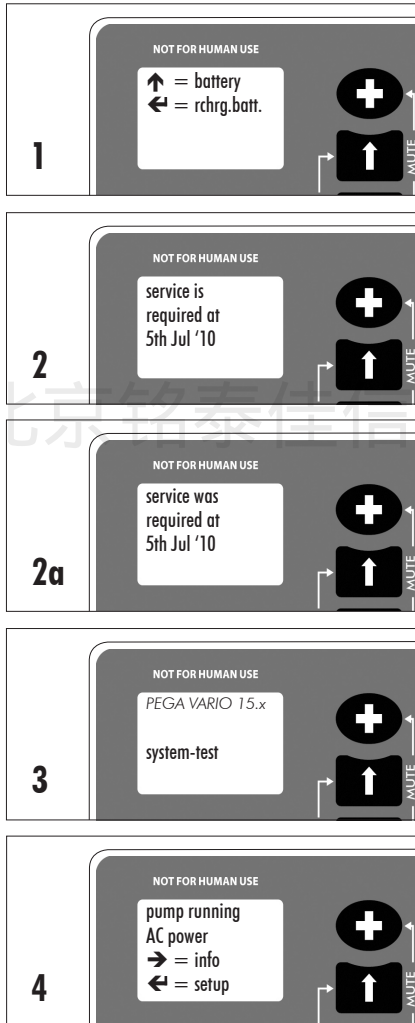
When power is supplied, ORCHESTA™ Model 500-V starts operating.

If no batteries or rechargeable batteries are fitted, the pump will stop as soon as the mains/AC plug is removed. Disconnect the mains/AC connection at the infusion pump first.

Remove the mains/AC plug

Take hold of the knurled metal ring (**see Figure 1**) and pull it back. The plug will be unlocked and can be easily removed.

Start-up / system-test



Start-up

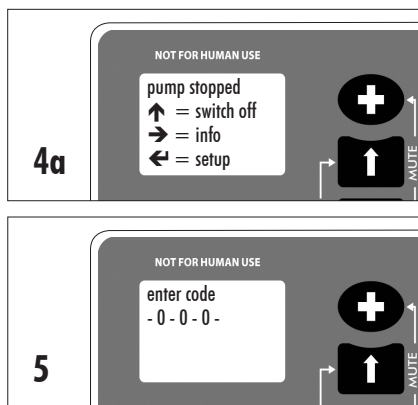
Procedure


1. If batteries or rechargeable batteries are inserted, **Figure 1** will automatically appear on the display.
2. Press the corresponding key. A short signal sounds.
3. A safety check must be carried out on the *ORCHESTA™ Model 500-V* at regular intervals (see p. 5). If the date is within the next 4 weeks, **Figure 2** will appear on the display.
4. If the date has already passed, **Figure 2a** appears on the display.

Caution

If the message appears or the date has already passed, send the pump immediately to the authorised technician or associate. If use of the infusion pump is continued beyond the next safety check date, the warranty expires and all product liability on part of Venner Medical (Deutschland) GmbH lapses.

5. **Figure 3** then appears on the display.
6. If any button is pressed during system test, "system-test interrupted" will be shown briefly on the display and **Figure 4** will appear, for example. The system-test will continue to run in the background.



7. If the pump stops and the power supply has been interrupted, the display of **Figure 4a** appears after reinstatement of the power supply and the test.
8. If a fault or warning occurs, **see Table 5, p. 46**.
9. If you want to change the programming, press the Enter/save  key. **Figure 5** appears on the display.
10. Enter the required code (see p. 25).

Additional system-test

The *ORCHESTA™* Model 500-V carries out a system-test on start-up (insertion of rechargeable batteries/batteries or connection of *PEGA®* POWER SUPPLY) and thereafter in the background at hourly intervals.

An additional system-test can be initiated by the staff.

Procedure

1. Disconnect the unit from the mains/AC.
2. Remove the batteries or the rechargeable batteries and re-insert them.

  **Warnings, messages**

Messages on the display

◆ **Rechargeable batteries or batteries weak:**

indication of weak batteries/
rechargeable batteries

◆ **Rechargeable batteries or batteries empty:**

red signal lamp lights up as soon
as delivery stops (**see p. 59**).

◆ **Type of power supply:**

additional information provided
in the “Pump running” state

北京铭泰佳信科技有限公司

Programming

All basic settings are made in settings menu 1 (Code 1 - 9 - 2). Settings menu 2 (Code 1 - 1 - 1) is only for setting the fluid supply when changing infusion accessories and when priming (for menu overview, **see p. 10-13**). The code protects against unintended or unauthorised changes to the settings.

ORCHESTA™ Model 500-V is supplied with the following default settings:

- ◆ Flow rate and supply: "0"

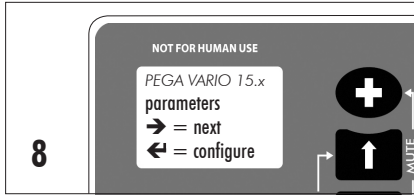
Calling up settings menu 1 - entering code



Procedure

1. Press the Numeric key once.
2. Press the Next key.
3. Press the Numeric key nine times.
4. Press the Next key.
5. Press the Numeric key twice (**Figure 7**).
6. Confirm the code by pressing the Enter/save key.
7. The message, "enter code - okay - ", will appear on the display.





8. After 3 seconds, “service required on...” will appear for a short time.

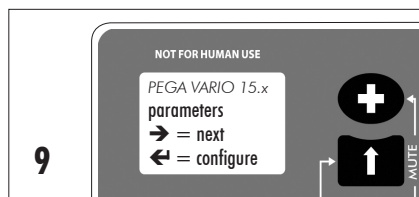
The pump is now in programming mode (**Figure 8**).

You can get to sub-menu B (Set) by pressing the Next key, and to sub-menu A (Configure) by pressing the Enter/save key.




北京铭泰佳信科技有限公司

Programming for continuous infusion






Procedure

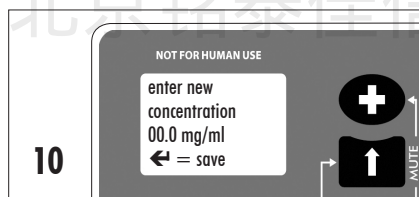
1. Call up code entry.
2. Enter code 1 - 9 - 2.
3. Select sub-menu A (Configure) **(Figure 9)**.





7. Repeated pressing of the Next key will lead you through the displays for setting the unit of measurement, the alarm signal, data transfer, pump clock and date. 

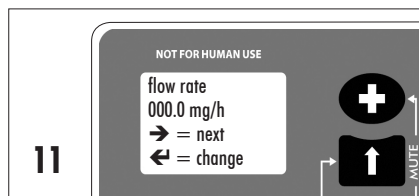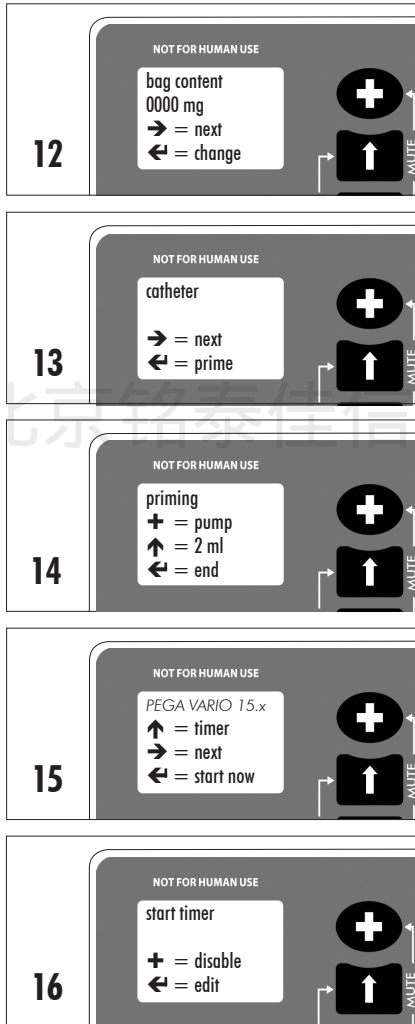
When you have completed sub-menu A (Configure) in settings menu 1, you automatically switch to sub-menu B (Set).
















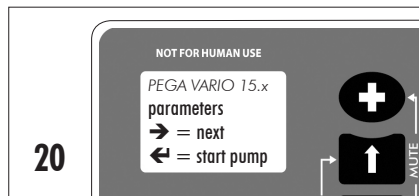
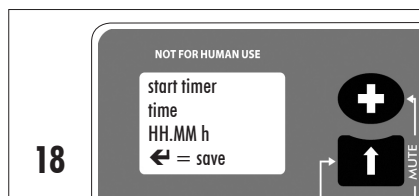
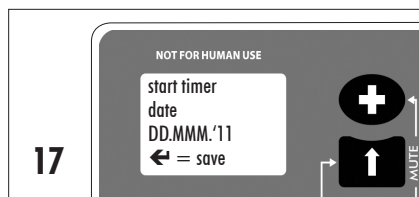
8. If you have set the pump for programming in mg/ml or $\mu\text{g/ml}$, enter the required value in the "concentration" window using the Numeric and Next keys **(Figure 10)**.  








9. Confirm the change with the Enter/save key to save the value. 
10. Press the Next key until you reach the menu item "flow rate" **(Figure 11)**. 
11. Set the required value using the Numeric and Next keys.  







12. Press the Enter/save key to save the new value. 
13. By pressing the Next key, you will reach the menu item “bag content” (Figure 12). 
14. Set the required value using the Numeric and Next keys. 

15. Press the Enter/save key to save the new value. 
16. By pressing the Next key, you will reach the menu item “catheter” (Figure 13). 
17. Press the Enter/save key in order to start the priming of the catheter system. 
18. Press the Special key until the catheter system has been primed or...
press the Numeric key. The pump will automatically deliver 2 ml (Figure 14). 

19. If you press the Enter/save key, Figure 15 appears. 
20. By pressing the Next key, you will reach the entry point of the menu (Figure 9); using the Up key, you will reach the menu item “start timer” (Figure 16). To start the pump press the Enter key. 








21. Press the Enter/save key in order to edit the date for scheduled pump start (**Figure 17**) or use the Special key to disable the timer.  

22. Set the date using the Numeric and Next keys and press the Enter/save key to save the new value. **Figure 18** appears.   

23. Set the time using the Numeric and Next keys and press the Enter/save key to save the new value. **Figure 19** appears.   


24. Set the required values for the stop timer in a similar manner. 

25. After timer values have been set, press the Next key until the menu item “parameters/start pump” appears (**Figure 20**). 

26. Press the Enter/save key to start the pump or the Next key to scroll through the settings again.  






27. If you want to delete parameters, proceed as described on **p. 29**.

Changing values

You can view the menu items by using the Next key. Settings can be changed individually or combined. To do so, follow the instructions on the display. 

Pay attention to the position of the cursor. The blinking cursor shows which numeral can be changed.

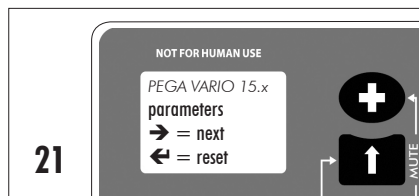
Procedure

1. Set the desired numeral using the Numeric key. 
2. Select the next numeral with the Next key. 
3. Set the desired numeral using the Numeric key. 
4. Repeat the procedure until all numerals have been selected.
5. Save your entries by pressing the Enter/save key. The saved value is displayed again. 
6. You can call up the next menu item by pressing the Next key. 

If date and time are changed, all parameters (if applicable, concentration, flow rate and supply) are automatically deleted.

Deleting all parameters


As the first step, any parameters set in the ORCHESTA™ Model




21



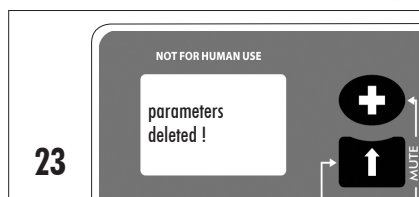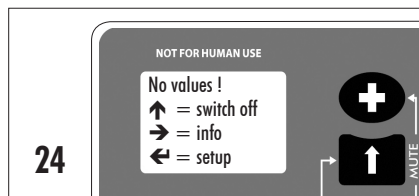22






23



24

500-V can be deleted. To do so, proceed as follows:


Procedure

1. Press the Next key. **Figure 21**  appears on the display.
2. Press the Enter/save key. **Figure 22**  appears on the display.
3. Confirm the following check by pressing the Numeric key ("yes"). **Figure 23**  Thus, all settings in settings menu 1, sub-menu A (Set), are reset to "0" (e.g. flow rate and if applicable, concentration). **Figure 23** appears.

Settings

As long as the pump is switched on, you can run through settings menu 1 as often as you like.

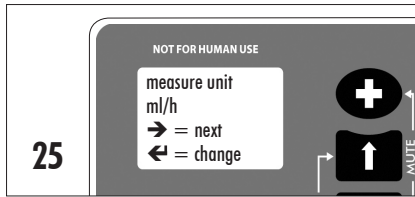
If "0" has been entered for flow rate, the pump identifies this to be an unreasonable entry. **Figure 24** appears on the display.


The ORCHESTA™ Model 500-V is switched off by pressing the Numeric key. 

Unit of measure

Procedure

1. Call up code entry.
2. Enter code 1 - 9 - 2.




3. Select sub-menu A (Configure).
4. You are now in the menu item "measure unit" (Figure 25). The display shows ml/h, mg/h, $\mu\text{g/h}$, mg/d or $\mu\text{g/d}$. The unit of measurement can be changed by pressing the Enter/save key. 




All further parameters are now shown in this unit of measurement. If you change between volumetric and gravimetric units, flow rate will be deleted.


北京铭泰佳信科技有限公司 Concentration

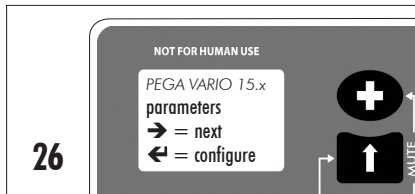If mg/h, mg/d, $\mu\text{g/h}$ or $\mu\text{g/d}$ has been selected as the unit of measurement, the concentration of the fluid must be entered. To do so, proceed as follows:


Procedure




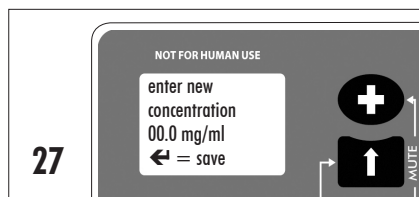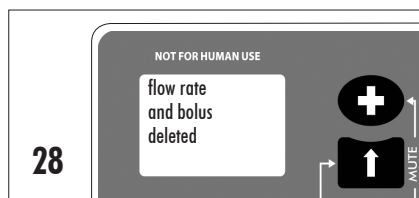1. Starting from the menu item "measure unit", press the Next key until **Figure 26** appears on the display. 
2. You can switch to sub-menu A (Set) by pressing the Next key. 
3. Press the Next key until you reach the menu item "Concentration". 




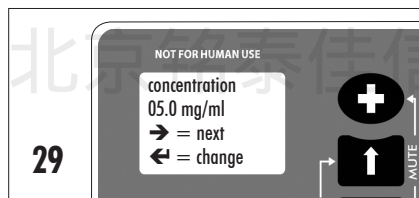
27

4. Enter the concentration using the Numeric and Next keys  (Figure 27).



28

5. Complete the entry by pressing the Enter/save key. **Figure 28**  appears briefly on the display.



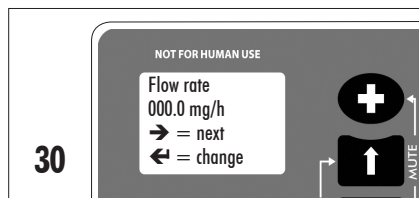
29

6. The new value is then displayed again, to check for correct entry. **Figure 29** appears on the display.






ORCHESTA™ Model 500-V automatically computes the corresponding volume.


Setting the flow rate

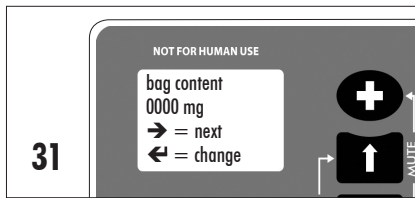Procedure







30

1. Call up code entry.
2. Enter code 1 - 9 - 2.
3. Select sub-menu B (Set).
4. Press the Next key until you reach the menu item "flow rate" **(Figure 30)**. 
5. Press the Enter/save key to enter the flow rate. 
6. Set the required value using the Numeric and Next keys.  
7. Press the Enter/save key to save the new value. 



Setting supply

Procedure

1. Call up code entry.
2. Enter code 1 - 9 - 2.
3. Select sub-menu B (Set).
4. By pressing the Next key, you will reach the menu item "bag content" (**Figure 31**). 
5. Set the appropriate value using the Numeric and Next keys.  
6. Press the Enter/save key to save the new value. 

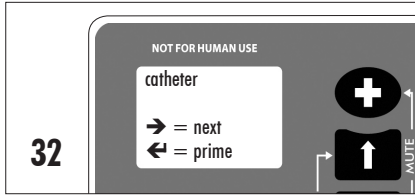
Priming

The infusion system must be primed before it is connected to the animal. The operation of the pump will be checked at the same time.


The priming function is available both in settings menu 1 (Code 1 - 9 - 2) and in settings menu 2 (Code 1 - 1 - 1).

Procedure

1. Call up code entry.
2. Enter the appropriate code.
3. If you are in settings menu 1, select sub-menu A (Set).



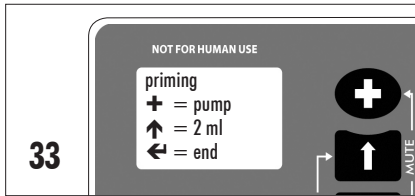4. Press the Next key in the settings menu until **Figure 32** appears on the display.




5. Press the Enter/save key. **Figure 33** appears on the display.



6. Press the Special key until the catheter system has been primed or...



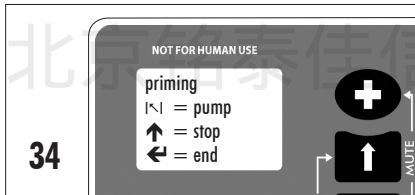
press the Numeric key. The pump will automatically deliver 2 ml (**Figure 34**).



7. Press the Enter/save key. **Figure 35** appears.



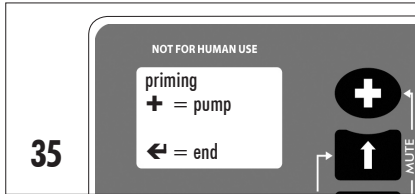If you press the Enter/save key again, it will take you back to **Figure 32**.




Leaving settings menu 1 - starting pump


Procedure

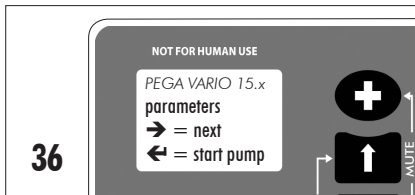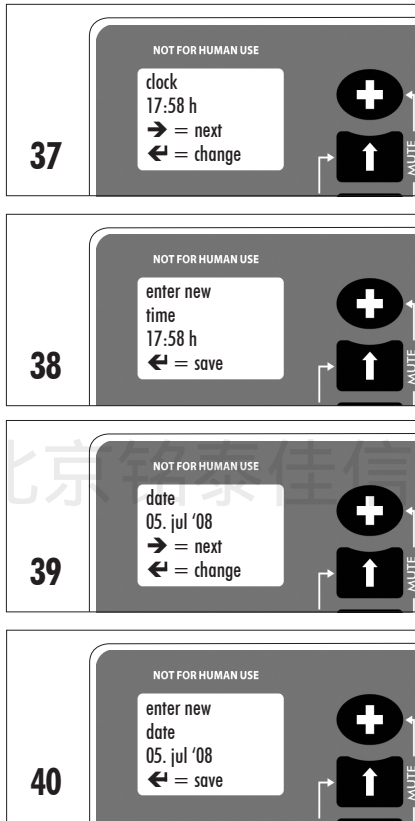1. Press the Next key until **Figure 36** appears on the display.












2. Press the Enter/save key to start the pump.





Setting pump clock and date

Procedure

1. Call up code entry.
2. Enter code 1 - 9 - 2.
3. Select sub-menu A (Configure).
4. Press the Next key until you reach the menu item “clock” **(Figure 37)**. 
5. Press the Enter/save key to change the time. 
6. Set the appropriate value using the Numeric and Next keys **(Figure 38)**. 

7. Press the Enter/save key to save the new time. 
8. Press the Enter/save key to change the date **(Figure 39)**. 
9. Set the appropriate value using the Numeric and Next keys **(Figure 40)**. 

10. Press the Enter/save key to save the new date. 

Caution

When the time and date are changed, the parameters in sub-menu B (Set) are deleted and must be re-entered. This applies to flow rate and concentration.

Brief instructions for daily care

Replace **PEGA® BAG** and **PEGA®** disposables

Fluids and disposables must be inserted when starting a new infusion or regularly replaced when changing and resuming an

existing infusion. The associated settings are made in settings menu 2 for technical staff (Code 1 - 1 - 1).

Warning

- ◆ It is essential to observe the general safety instructions (**see p. 5**).
- ◆ Check that the pump is working and the unit is in proper condition each time before using it.

To do this, a self-test must be carried out when starting up (**see p. 22**) and the pump must be primed (**see p. 44**).

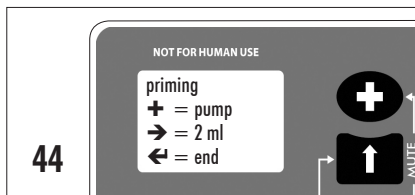
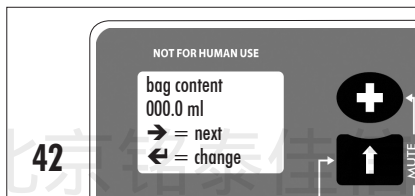
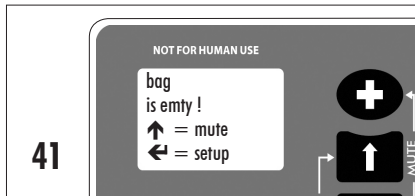
- ◆ To avoid contamination, carry out all operations using proper aseptic techniques. Disinfect your hands before starting work and avoid touching the open Luer connections.
- ◆ Only use fluids or solutions which are compatible with the material of the disposable items.

Note

- ◆ If the fluid supply does not exceed two hours, the time remaining until the contents of the bag are exhausted is shown on the display.









Between 8 a.m. and 8 p.m. this message is signalled every hour with a beeping tone (this is an artifact from human use).

- ◆ When the pre-set supply has been used up, the red signal lamp of the **ORCHESTA™** Model 500-V lights up, an alarm signal sounds and **Figure 36** appears on the display. By stopping the pump this alarm can be cut off and the pump may be switched off until it is used again (**see p. 40**).



Setting supply and priming

Procedure

1. Connect a pre-filled PEGA® BAG.
2. Press the Enter/save key. 
3. Enter code 1 - 1 - 1, **see p. 12.**
4. Press the Enter/save key.  **Figure 42** appears on the display.
5. Press the Enter/save key.  **Figure 43** appears on the display.
6. Set the supply using the Numeric and Next keys.  
7. Confirm your entry by pressing the Enter/save key. 
8. Press the Next key and **Figure 44** will appear on the display. Prime the infusion system. Exit settings menu 2 (Code 1 - 1 - 1) by pressing the Enter/save key.  

i Note

- ◆ There is a restriction to 5 ml for 10 hours when priming in settings menu 2 (Code 1 - 1 - 1). If this volume is reached, the priming function is blocked.
 - ◆ The priming process only works if a supply is set before priming.
 - ◆ Care must be taken to ensure that there are no air bubbles in the catheter and the filter.
- ◆ If air bubbles are still seen during the priming process, gently flick your finger against the filter and the tube until all air bubbles have been flushed out.
 - ◆ Continue the priming process until all air bubbles have been removed from the line.

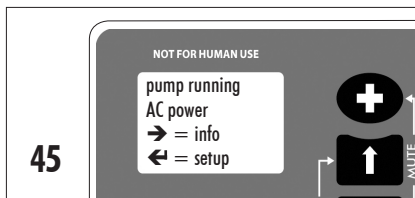
北京铭泰佳信科技有限公司

Connection with animal and starting the infusion

If all the settings have been entered and preparations for the infusion completed, *ORCHESTA™* Model 500-V can be connected to the animal and started.

Procedure

1. Connect the catheter system and pump to the animal.
2. If you have called up settings menu 2 (Code 1 - 1 - 1) when the pump was running, the pump will start as soon as priming is completed (**Figure 45**).
If you have called up settings menu 2 (Code 1 - 1 - 1) when the pump was stopped, you will have to start the pump. Exit the menu by pressing the Enter/save key. **Figure 46** appears on the display. To start, **see p. 40**.
3. **Figure 47** will appear if the set supply is sufficient for 2 hours or less.




After starting, the *ORCHESTA™* Model 500-V pumps at the set flow rate.



The green signal lamp blinks at regular intervals.

Pressing any key will take you back to the information menu.

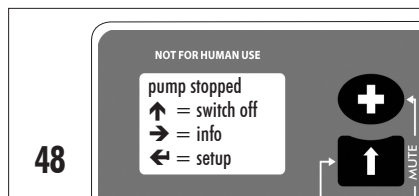
Stopping / switching off / re-starting the pump

Stopping and switching off the pump

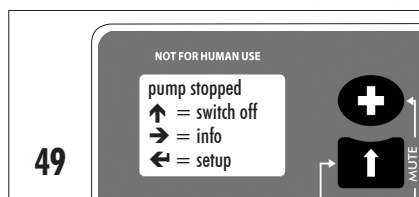


1. Press the Numeric key and the Enter/save key simultaneously.   **Figure 48** appears on the display. The pumping process will stop.




2. Press the Numeric key again.  **ORCHESTA™ Model 500-V** will now switch off (**Figure 48a**). The values set remain saved.


Re-starting the pump



1. Press any key. **Figure 49** appears on the display.



To check the set values, press the Next key. The information menu will appear (**see p. 42**). 



2. Press the Enter/save key to enter the code 1 - 1 - 1: **Figure 50**. 

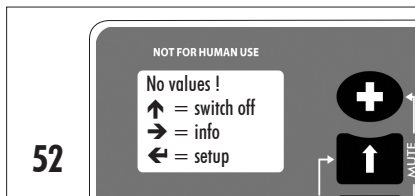
If you do not press a key, the screen will switch off again after a short time.



3. Press the Special key. **Figure 51** appears on the display. 
4. Press the Enter/save key and the pump will start. 

Switching off the pump when infusion parameters have not been programmed

If you want to switch off the pump for an extended period without entered values, delete all parameters in Settings menu 1 (Code 1 - 9 - 2) (**see p. 29**).



Note

If mg/h (mg/d, μ g/h or μ g/d) has been set for the pump in sub-menu A, a concentration must be entered before switching off.

Figure 52 will appear on the display when restarting.

Information menu

When the pump is operating all current values can be called up and checked at any time in the information menu. The information menu differs depending on whether the unit of measurement has

been set to ml/h or mg/h or $\mu\text{g/h}$, or mg/d or $\mu\text{g/d}$.

By pressing the Next key, you will reach the following items in succession: **see Table 2.**



	ml/h - version	mg/h, $\mu\text{g/h}$, mg/d or $\mu\text{g/d}$ - version
1a	Supply/flow rate	Concentration/flow rate
1b		Concentration/supply
2	Range of supply	Range of supply
3	Amount of fluid administered since date/time	Amount of fluid administered since date/ time
4	Start timer	Start timer
5	Stop timer	Stop timer

Table 2: Information menu

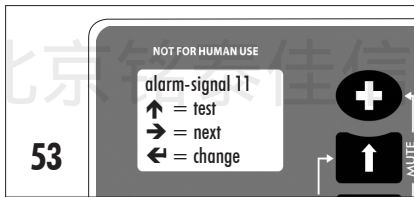
Alarm and warning functions

Setting alarm signals

Users can choose different tonal sequences to ensure that audible alarm signals can be heard more easily in the daily routine.

Procedure

1. Call up code entry.
2. Enter code 1 - 9 - 2.
3. Select sub-menu A (Configure).
4. Press the Next key until the "Set alarm signal" menu item appears (**Figure 53**).
5. Press the Enter/save key until you reach the desired tonal sequence.
6. Press the Numeric key to test the selected tonal sequence.
7. Press the Next key to confirm the selected signal tone.



Note

Internal warning tones cannot be changed.

Signals

Audible signals

- ◆ Every press of a key is followed by a beeping tone which confirms that the action was successful.

In the event of malfunctions, signals sound in conjunction with a corresponding warning message on the display.

- ◆ In settings menus 1 and 2 a short beeping tone sounds every 15 seconds.

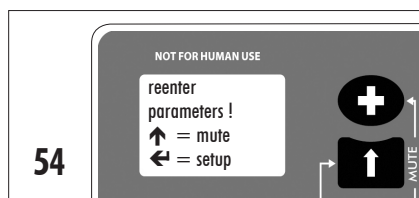
Audible alarms

北京铭泰佳信科技有限公司 You can find the meaning of the audible alarms in the table below.

The audible warning can be suppressed for approx. 5 min in each case.

Procedure

1. For example, **Figure 54** appears on the display.
2. Press the Numeric key. The audible alarm is suppressed for approx. 5 min (the pump will not work during this time!).



Audible alarm	Significance
Intermittent signal with 1 s tone / 1 s pause	No fluid supply left
Intermittent signal with 0.5 s tone / 1 s pause	Occlusion
Intermittent signal with 1 s tone / 10 s pause	All other alarms

Table 3: Audible signals

Mute/unmute audible alarms

For muting/unmuting all key tones and audible alarms, press Next key and Enter/save key in the Information menu at the same time. This setting will remain unchanged in case of switching off the device or changing the batteries.



Audible alarm when power supply is interrupted (if not muted)

If the power supply is interrupted while the pump is running, an alarm tone (with lower frequency) sounds which differs from the above mentioned tones. This alarm can continue for up to 15 min if the fault is not corrected.

Visual signals

Audible and visual signals are complemented by corresponding messages and symbols on the display (exception: if the power supply is interrupted, there is only an audible alarm).

Visual signal	Significance
Signal lamp flashes green	Pump running as programmed
Signal lamp not flashing	Pump not running
Signal lamp shining red	There is a fault and the pump is not running

Table 4: Visual signals

Tracing and correcting faults

ORCHESTA™ Model 500-V offers maximum safety as the pump has built-in checking functions which constantly monitor the infusion settings. As soon as a malfunction occurs, the pump will set off a visual and audible alarm.

The following table shows which messages can occur and how you can correct the malfunctions.

⊗ Red signal lamp

i Note

If you are unable to rectify a fault, send the pump and if necessary the PEGA® POWER SUPPLY in the pump case together with the fitted consumables to the authorised dealer. When dispatching, please observe the regulations for packing potentially contaminated goods.

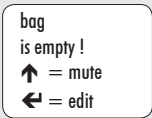
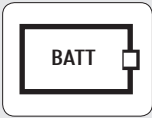
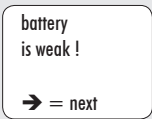
Display	Cause	Correction
1.  ⊗	PEGA® BAG empty	Change infusion system and restart pump or stop pump and switch off
2. 	Limited life in batteries	Insert new batteries or new rechargeable batteries or connect pump to mains/AC
3. 	1. Limited life in batteries on start-up (max. 20 % of capacity left) 2. Rechargeable batteries declared as batteries by mistake	1. Insert new batteries 2. Refit rechargeable batteries and declare them correctly; connect mains/AC if necessary

Table 5: Fault list

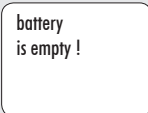
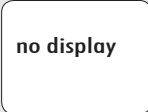
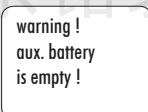
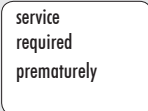
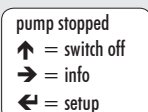
Display	Cause	Correction
4. 	Batteries used up. Infusion can stop any minute; not possible to start new infusion	Insert new batteries or new rechargeable battery pack or connect pump to mains/AC
5. 	1. Pump is switched off 2. No batteries or batteries empty; alarm tone (with lower frequency) for up to 15 minutes if power interrupted while pump is running 3. Technical fault in pump	Activate display and pump with any key or insert new batteries or rechargeable battery pack or connect pump to mains/AC; if this does not rectify the fault, return unit to authorised dealer
6. 	No internal storage power	Send pump back to authorised dealer
7. 	Service required for technical reasons (e.g., max. delivery volume reached)	Send pump back to authorised dealer
8. 	Pump has stopped	Check values and re-enter if necessary; restart or switch off pump

Table 5: Fault list

Display	Cause	Correction
9. <div style="border: 1px solid black; border-radius: 10px; padding: 5px; width: fit-content;"> reenter parameters ! ↑ = mute ← = edit </div>	When working in settings menu 1: - No key pressed for 30 min - The batteries have been removed	Check values in settings menu 1 and if necessary, re-enter
10. <div style="border: 1px solid black; border-radius: 10px; padding: 5px; width: fit-content;"> system error ! ↑ = mute ← = edit </div>	Pump has detected corrupted data in internal test	Call up settings menu, check settings and re-enter values
11. <div style="border: 1px solid black; border-radius: 10px; padding: 5px; width: fit-content;"> catheter occlusion ! ↑ = mute ← = edit </div>	Occlusion on animal's side (blockage): Needle / filter is blocked or catheter system kinked	Check catheter system and remove occlusion, change catheter system if necessary
12. <div style="border: 1px solid black; border-radius: 10px; padding: 5px; width: fit-content;"> please check pump head ! ↑ = mute ← = edit </div>	Check PEGA® PUMP HEAD V; fault in pump/PEGA® PUMP HEAD V	Change PEGA® PUMP HEAD V
13. <div style="border: 1px solid black; border-radius: 10px; padding: 5px; width: fit-content;"> please close pump ! ↑ = mute ← = edit </div>	Housing clips have not clicked, home or no PEGA® PUMP HEAD V fitted	Fit PEGA® PUMP HEAD V, close lid, engage pump lid clips, re-start pump

Table 5: Fault list



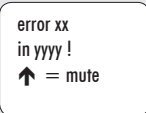
Display	Cause	Correction
14. 	Internal pump fault (coded fault message) Technical fault in pump	Switch off power supply, re-start pump; if fault is not corrected, send pump to the authorised dealer
15. 	Technical fault in pump	Switch off power supply, re-start pump; if fault is not corrected, send pump to the authorised dealer

Table 5: Fault list

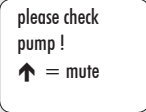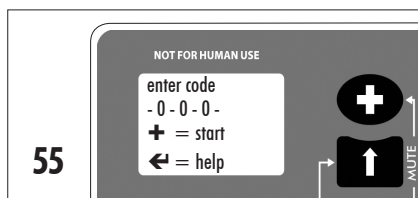Re-starting the pump






Procedure

1. If possible, correct the cause of the alarm.

2. Press the Enter/save key. 

Figure 55 appears.

Pressing the Enter/save key again causes a help text to be displayed. 

3. Press the Special key and then the Enter/save key. The pump is re-started.  

北京铭泰佳信科技有限公司

What should I do if *ORCHESTA™* Model 500-V is dropped?

Examine the pump, pump lid and *PEGA®* BOX for external damage (cracks, etc.).



If the pump is damaged, send it to the authorised dealer.

Also check the disposables for damage and replace them if necessary.

If the *PEGA®* BOX is damaged, replace it.

If the pump is outwardly undamaged, proceed as follows:

Procedure

1. Press the Numeric key and the Enter/save key simultaneously (Stop).  

2. Briefly interrupt the power supply (pull out the plug or remove batteries/rechargeable batteries and then re-insert them).

ORCHESTA™ Model 500-V will now carry out a system-test. If the system-test is successful, the pump is ready to operate and can be re-started. Otherwise, send the pump to the authorised dealer.

Cleaning and safety check

Cleaning

- ◆ Wipe the *ORCHESTA™* Model 500-V with a damp cloth.
- The pump is resistant to alcoholic disinfectants.

Safety check

- ◆ For the *ORCHESTA™* Model 500-V, this safety check must be carried out by a trained technician at least every five years or within one year after a delivery volume of 50 liters.
- ◆ The date of the next safety check is displayed in the settings menu. Within the final 4 weeks before the next service, this is notified on the display every time the batteries are inserted.
- ◆ If the service is due or the date has been exceeded, send the *ORCHESTA™* Model 500-V immediately to the authorised dealer.
- ◆ If use of the pump is continued beyond the safety check date, the warranty expires and all product liability on the part of Venner Medical (Deutschland) GmbH lapses.

Dosing accuracy confirmation

Basal rate	Time
15 ml/h	24 s
10 ml/h	34 s
5 ml/h	64 s
1 ml/h	6 min
0.1 ml/h	60 min

Tab. 6: Delays before occlusion alarm

Dosing accuracy is continually checked and registered by a dual monitoring system of two microprocessors within the pump.

This ensures that an alarm is triggered at an inaccuracy of 75 μl (**check PEGA® PUMP HEAD V**).

Table 6 shows the maximum delay times before the occlusion alarm is activated.

An occlusion on the animal's side results in a bolus volume in the infusion system. After the occlusion has been removed, this volume is dispensed. The values in **Table 7** have been determined for the following PEGA® disposables:

Designation	Max. bolus volume if occlusion 1,800 hPa
PEGA® LINE 100 SF + PEGA® PUMP HEAD V	60 μl

Tab. 7: Bolus volume for occlusion

The infusion systems are designed for a pressure $p \leq 2,000$ hPa (29 psi).

Caution

The use of consumables not approved by Venner Medical (Deutschland) GmbH, can endanger the animal, e.g.:

- no additional security against free-flow in the event of a malfunction

- no protection against air infusion
- larger bolus volume
- overdosage or underdosage of medication without alarm
- pressure stability not sufficient
- leaking of disposable items

Compliance

This manual has been written with consideration to the requirements in the International Standard IEC 60601-1 for Medical Electrical Equipment.

The *ORCHESTA™* Model 500-V infusion pump is designed to comply with (among others) the following international safety and EMC requirements for medical electrical equipment:

- ◆ IEC 60601-1 (1988 with Amendments) medical electrical equipment, which classifies the *ORCHESTA™* Model 500-V as follows:
see table 9, page 54.
- ◆ IEC 60601-1-2 (2001): Electromagnetic compatibility

北京铭泰佳信科技有限公司

Technical data and accessories

Technical data

Classification	Type BF
Housing	ABS, shock proof, splash proof protection degree IPX 4
Dimensions	87.0 x 64.2 x 33.0 mm
Weight	152 g, without rechargeable battery pack
Drive	Micro double-piston system
PEGA® RECHARGEABLE (rechargeable battery pack)	2.4 V, NiMH
Batteries	2 x 1.5 V, Mignon (R6, AA)
Storage temperature	-30 °C bis +40 °C
Operating temperature	+5 °C bis +40 °C
Rel. humidity	between 30 % and 75 %
Flow rate	0.01 ml/h to 15 ml/h adjustable in steps of 0.01 ml/h, similar setting options when programming in mg/h or µg/h
measure unit	ml/h, mg/h, µg/h, mg/d, µg/d
Supply	0.1 ml to 999.9 ml adjustable in steps of 0.1 ml, similar settings in mg and µg
Accuracy	± 2 % or ± 2 µl/h
Pressure sensor	Cut-out pressure for occlusion 1,600 ± 200 hPa (23 ± 3 psi)
Max. delivery volume until next service check due	50 l, if last service check more than one year ago
Data transmission	Compatible with ORCHESTA™ automated infusion system

Tab. 8: Technical data

Delivery contents *ORCHESTA™* Model 500-V

Number	Contents
1	<i>ORCHESTA™</i> Model 500-V
2	Batteries (1.5 V, Mignon, R6 AA)
1	Pump case

Table 9: Delivery contents

Further accessories, disposables

Accessories/disposables are listed on a separate leaflet.

北京铭泰佳信科技有限公司

Appendix

Warranty, claims and liability

Venner Medical (Deutschland) GmbH offers warranty against defects and damages as part of the following rules:

The warranty period is 24 months from delivery of the infusion pump by Venner Medical (Deutschland) GmbH. If there is a defect, Venner Medical (Deutschland) GmbH is entitled to choose between correcting the fault or replacing the unit.

If Venner Medical (Deutschland) GmbH is not prepared or able to correct the fault/replace the unit, in particular if the process of doing so exceeds reasonable time limits for reasons for which Venner Medical (Deutschland) GmbH is responsible, or if the attempt to correct the fault has failed twice, the buyer shall be entitled to choose between withdrawing from the contract or demanding a corresponding reduction in the purchase price.

Claims for compensation arising from breach of contract or unauthorised action are not permitted either against Venner Medical (Deutschland) GmbH or their various agents unless there is a case of deliberate or grossly negligent conduct. This also applies to claims for damages due to non-performance, but only in so far as compensation is demanded

for indirect or consequential damages unless liability is based on the warranty which is intended to protect the buyer from the risk of such damages.

In the event that an important contractual obligation is not met, liability shall be limited to the foreseeable damages when the contract was concluded.

In order to claim the right to have a defect corrected or a replacement unit supplied, it is necessary to send the defective infusion pump to Venner Medical (Deutschland) GmbH or the responsible dealer. The buyer shall bear the carriage costs.

The warranty obligation shall not extend to faults which are due to...

- ◆ extraordinary or unavoidable events
- ◆ repairs carried out by unauthorised persons or to third party intervention
- ◆ improper operation and handling, in particular improper storage, dropping, immersion in liquid, unintended use and any disregard of the instructions for use
- ◆ the use of disposables that are not approved by Venner Medical (Deutschland) GmbH

Trumpet curve (accuracy diagrams)

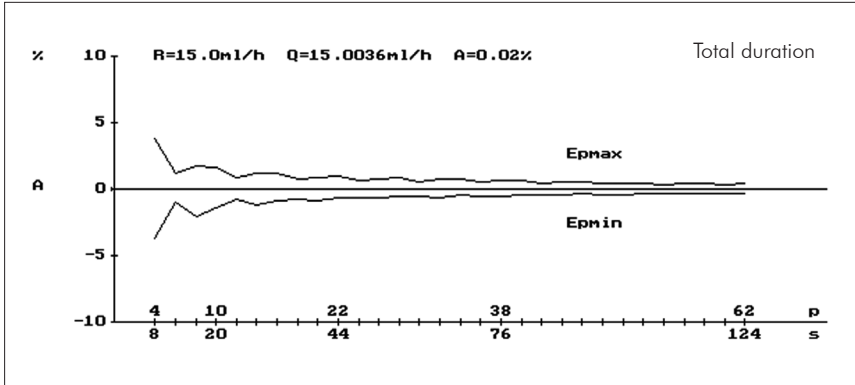


Figure J: Trumpet curve

The coordinate system of the trumpet curve shows the percentage deviation from the set flow rate on the y axis. The x axis shows the observation period. The possible deviation from the set flow rate can be seen from the trumpet curve for any observation period.

For example, this means for a flow rate of 15 ml/h:

in an interval of 2 minutes, the maximum deviation, E_{pmax} for the PEGA® vario is approx. 0.5 %, and the minimum deviation, E_{pmin} approx. 0.5 %. The total error is 0.0 %.

Index

A

accessories 4, 5, 55
additional system-test 22
air bubbles 38
alarms
 audible 44
 visual signals 45
alarm signals
 setting up 43

B

background lighting 9
batteries 23
battery compartment
 latch 8
 lid 19
beeping tone 36, 44
bolus functions 10
brief instructions for daily care 36

C

case 14
catheter system 34, 39
claims 56
codes 10, 11, 12
 entering 24
connection 20
contamination 36
continuous infusion 26
correcting faults 46

D

date
 deleted parameters 29
 setting 35
defects and damages 50, 56
delay times 52
deleting
 all parameters 41

disinfectants 51
display 8
disposables 4, 36, 55
drive pinion 16

E

enter/save key 14

F

fitting aid 16
flow rate 30, 32, 39, 42
fluid supply 36
function keys 8

G

gravimetric units 31

H

hardware 8

I

information menu 9, 42
infusion
 starting 39
infusion system 37

L

liability 51, 56
luer-lock connection 4

M

magnetic fields 6
mute/unmute 14, 45

N

next key 14
numeric key 14

O

occlusion 52, 54

P

parameters 29

- delete 41
- PEGA® BAG 4, 5
- PEGA® BOX 8
- PEGA® BOX and soft cases 17
- PEGA® catheter systems 4, 5
- PEGA® disposables
 - replacing 36
- PEGA® POWER SUPPLY 9
- PEGA® PUMP HEAD V
 - inserting 16
- PEGA® reservoir 4
 - replacing 36
- polarity 19
- power supply 8, 23, 45
- power supply unit 18, 20
- pressure sensor 17
- priming 33, 37
- product liability 21, 51
- programming 22
- pump
 - open 15
 - re-starting 40, 50
 - starting 36
 - stopping 40
 - switching off 40, 41
- pump lid clips 8, 15
- R**
 - rechargeable battery pack/batteries 18
 - inserting 19
- S**
 - safety check (SC) 51
 - safety Information 5
 - scope of delivery 55
 - security check 30
 - settings 29
 - settings menu 1 9, 11, 30
 - settings menu 2 12, 38
 - signal lamp/infrared interface 8
 - signals
 - visual 45
 - special key 14
 - start timer 11
 - stop timer 11, 28
 - storage 14
 - sub-menu A/B 10
 - supply 33, 36, 37, 39
 - switching off 40
 - system-test 21, 22, 50

T

 - technical data 54
 - time
 - deleted parameters 29
 - tonal sequences 43
 - tracing and correcting faults 45
 - transport and storage 14
 - trumpet curve 57

U

 - unit of measure 13
 - unit of measurement 10, 12

V

 - visual signals 45
 - volumetric units 31

W

 - warning tones 43
 - warranty 21, 51, 56

X

 - x-rays 6

Z

 - zinc carbon batteries 6, 18

Meaning of the symbols on the rating label and on the packaging

	Expiry date
	Item for single use
	Caution, consult accompanying documents
	Batch
	Serial number
	Date of manufacture
	Sterilisation with ethylene oxide
	Sterilisation by radiation
	Compliant with product-specific European Directives
	Purchase code number
	Open here
	Protect against humidity
	Storage temperature
	Classification BF
	Recycling
	Dispose of pump, rechargeable batteries and batteries in accordance with the statutory regulations
	Protect against direct solar irradiation

Concept & Realisation: www.arno65.de

Mar 2011

Open the pump lid to see the rating label (**see p. 15**)

manufactured by:

exclusive lab animal distribution by:

VENNER
Medical (Deutschland) GmbH

Mühlenstraße 17
D - 24229 Dänischenhagen b. Kiel
Phone +49 (0) 4349 91 54 0
Fax +49 (0) 4349 91 54 10
info@vennermedical.de

INSTECH

Instech Laboratories, Inc.
5209 Militia Hill Road
Plymouth Meeting, PA 19462 USA
www.instechlabs.com
Tel (800) 443-4227
Tel +1-610-941-0132
Email support@instechlabs.com

ORCHESTRA™ Model 500-V
E 4002001A