

CMA 120 System for Freely Moving Animals

INSTRUMENTS



CMA 120 System for Freely Moving Animals

CMA 120 Plastic Bowl with Food and Water Containers

- 带有双通道转环的平衡比
- 确保管道远离动物，防止扭曲



订购信息

CMA 120 System for Freely Moving Animals

Ref. No.

CMA 8309049

Includes:

- CMA 120 Plastic Bowl
- CMA 120 Swivel Assembly
- CMA 120 Balance arm
- Plastic Collar 100 pcs
- Tubing Adapters 10 pcs
- FEP tubing 1 m
- Vial, Plastic 300 μ L, 25 pcs
- Caps, Plastic, 25 pcs
- CMA 120 Plastic Bowl with Food and Water Containers

CMA 8309672

The CMA 120 System for Freely Moving Animals

支持清醒状态下长时间的小型实验动物的微透析研究. CMA 120仪器可与任何一种微透析系统结合使用. 微透析探针连接至CMA注射泵、CMA 110液体开关, 并通过双通道旋转接头连接至任何CMA 142或CMA 470收集设备. 转环安装在平衡臂上, 允许动物自由移动. 旋转支架用于固定一根带有套环连接器的导线和两个用于300 μ L收集瓶的固定器. 连接到动物项圈的钢丝转动旋转接头并支撑管子.


当植入两个微透析探针时, 或当微透析与通过旋转接头的一个通道局部注射相结合时, 使用手动操作收集.

The CMA 120 Bowl with Food & Water Containers

用于将自由移动的动物放置更长时间的研究. 添加这些容器后, 动物可以随意进食和饮水. 食物容器和水瓶设置在碗的外侧, 以避免在碗内时干扰拴系动物的运动. 这两个容器都可以很容易地拆下进行清洁和重新填充.

CMA 120 System for Freely Moving Animals

Accessories

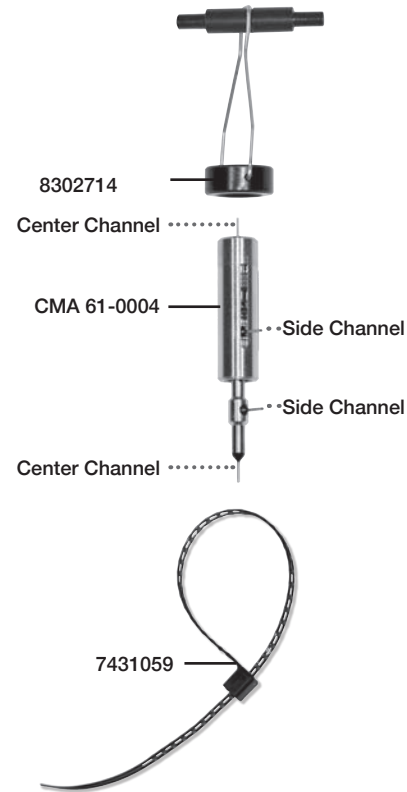

INSTRUMENTS

ORDERING INFORMATION

CMA 120 System for Freely Moving Animals, Accessories

	Ref. No.
CMA 120 System without Bowl	CMA 8409029
CMA 120 Plastic Bowl	CMA 8309031
CMA 120 Swivel Assembly	CMA 8309048
CMA 120 Swivel Assembly without Swivel	CMA 2409090
CMA 120 Wire Set with Collar Connector	CMA 2409051
CMA 120 Balance arm	CMA 8309032
Gimbal for Swivel	CMA 8002714
Dual Channel Swivel	61-0004
Dual Channel Swivel for mouse	72-0000
Clamp for In Vitro Holder	CMA 8409033
Plastic Collar 100/pkg	CMA 7431059
Tubing Adaptors 10/pkg	CMA 3409500
FEP tubing, 1 m	CMA 3409501
FEP tubing, 10 x 1m/pkg	CMA 8409501

CMA 120 Instrument Table

For use in organizing all the components of your micro-dialysis experiment for ease of operation

CMA 8309046

See accessories on next page.